

Not your average
manual



burstTM
WATER FLOSSER

Hamish & Brittany

Best friends & founders of BURST

Now streaming

Hey, fellow water flosser!

It's been a long time coming — thousands of hours of research, development, and more professional testing than they do on an airplane (well, probably). All hail, water flosser royalty has landed!

Some people think water flossers are more funky than fresh. No doubt, they've earned a funny rep. Many are the size of a small chair, look like a robot that might turn on you at any minute, and have more wires than a utility pole. Not the BURST Water Flosser! It's as modern, compact, and stylish as they come. It's also cordless!

Got tightly spaced teeth? No worries! With its unbelievable power and precision, our flosser easily cleans those hard-to-reach areas of the mouth. Interdental spaces and periodontal pockets — places where bacteria hang out — are a cinch to get into. And with your mouth free of food and debris, you'll be fresh all day.

All in all, it's super easy preventative oral care. This one-way ticket to Fresh Breath City has been developed with and offered by thousands of BURST Ambassador Dental Professionals. Your dentist will love you the next time they see you. We're sure they won't miss using dreaded words like root canal, cavity, or gum disease!

Brittany & Hamish

Best friends & founders of BURST



Ladies and gentlemen, this is your captain speaking

On behalf of the flight crew, I'd like to welcome you aboard BURST Flight WF to Shiny Town. This is Captain 'Go with the Flow' (yes, I am punny) and I'll be taking you through our user manual today!

1

Getting started with your BURST Water Flosser

First up, click your chosen tip into the top of the neck. You can rotate it 360 degrees to get at any angle you like.

2

Fill me up, buttercup!

Don't forget, this plane needs fuel. In this case, water! Lift the silicone cover on the back of the reservoir, and fill it up with room temperature or warm water — but not hot! Then, position the tip along your gumline and press the happy smiley face. It's showtime!



As told by Wing Commander,
Dental Hygienist Krismichelle Waller

The technique AKA how to fly the plane!



Step 1


No splash zone

Make sure you're over the sink (or better yet, in the shower) to spare your bathroom floor from a surprise splashing!

Step 2

Extra love here

Aim the jet tip at the gumline at a 90-degree angle, then press the smiley to start the engine. Move the tip along your gumline, giving extra love to the spaces between your teeth.

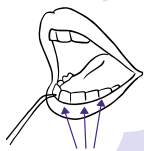

Extra love here

Step 3

Sparkle



Sparkle!

After about 20-40 seconds, refill the reservoir and keep on flossin' until you've completed your entire mouth.

Tip one

Your flosser is water-resistant so don't be afraid to use it in the shower!

Tip two

Keep it upright so you can be sure to thoroughly deplete the water reservoir.



Tip three

Refill as many times as needed to get that whole mouth clean.

Tip four

Need a break?
Mouth filling with water?!
Press the smiley again to pause the jet and catch your breath.

Pro tips



Krismichelle Waller
Dental Hygienist
Wing Commander

Changing modes

We're not ones to boast, but in this instance, it's the absolute truth: our engine is VERY powerful. That means our default mode is **standard** (the middle setting).

Once you've racked up some decent flossing hours, you'll be ready for the highest setting — **turbo**!

Hold down the smiley face to change modes (handy lights will flash to indicate mode choice). Toggle from the top mode (**turbo**) to the middle setting (**standard**) and on to the bottom (**pulse**).

● Turbo

Full speed — like Mach 3, only better!

○ Standard

The place to start; totally powerful for a squeaky clean.

⊗ Pulse

Bit of a massage mode; a great tidy-up with a little less oomph.

+ Red light, don't panic!

This is not a drill but be cool — we're only low on charge! Just insert the charging cable in the teeny-weeny hole halfway up the back (above the water reservoir). Then sit back and enjoy the inflight entertainment — we'll be fully charged in six hours.



Switching tips is easy-peasy



No need to call the flight attendant, just push the eject button on the back of the neck (NOT under your seat) and remove the tip.

Your BURST Water Flosser comes with the classic tip for daily use. There are two other specialist tips available for specific needs. Just go to www.burstoralcare.com to nab them.



Classic tip

The everyday tip, the OG, Old Faithful, the full-mouth plaque blaster. We'll call it the all-arounder. You can't go wrong with this one. Unless you have a specific condition or dental work, this is the tip to rely on.



Perio tip

Gentle is the name of the game here. A slim nozzle with a slight bend delivers water into crevices and problem spots. Place the perio tip at a 45-degree angle in the desired area and move as needed. You can always toggle on/off to hit particular pockets and save fuel!



Ortho tip

This tip is built to clean in and around orthodontic wires and brackets. Its tiny clustered bristles navigate hard-to-reach places where metal gets in the way.

Note: Ortho is also great for brushing away plaque and food debris for you non-ortho folk.

Shout out to the BURST Ambassador Dental Professionals!

Our Pros are involved in so much of what we do, no foolies! Our water flosser is no different. From the word go, they work like mad; advising, designing, testing, tweaking, anything to make sure this bad boy is the best of the best! It took over a year of refinement (and a tear or two) over in BURSTlabs.

Our techies and Pros have worked their cotton blend socks off to create something powerful, but compact. The BURST Water Flosser is sleek and stylish without compromising efficiency. We think they've nailed it.

30k

**BURST Ambassador
Dental Professionals**



Refreshment plan

The super duper BURST Ambassador Dental Professionals recommend replacing your flosser tip every 24 weeks. That's because they want you looking and feeling fresh! As long as you're subscribed to your replacement tips, we'll proudly offer a lifetime warranty on your device and replace it free of charge. Woohoo!



Questions? We've got the answers

What if I lose my tip?!

Take it easy like Sunday morning. Simply log into your account and order a new one ahead of your next scheduled delivery. We'll let you know when it's coming!

Can I water floss in the shower?

Jumpin' Jiminy, of course you can! We also love long bubble baths and late-night strolls on the beach. We'd go anywhere with you...

How do I clean my water flosser and how often should I?

Unplug before you get down to business. You can remove the reservoir for easy cleaning using warm water (95 degrees — Fahrenheit, not Celsius!), and a cloth. Simply slide the tab at the bottom of the device to unlock the tank. Then, give the tank a little downward tug — it should come right off.

For the tip and the unit, use a dry cloth to rub off any grubby grub. Clean as often as needed.

Where and how do I store my device?

Somewhere clean and dry like your countertop or medicine cabinet. Make sure the reservoir is emptied after and between uses for storage.

How do I reset my water flosser?

All you gotta do is...

- Plug the water flosser into a wall charging block using the included USB cable.
- Hold the power button for three seconds.
- The + sign will flash red to show it has been reset.

Also, it could be that your flosser is out of juice! In that case, insert the charging cable into the teeny-tiny hole halfway up the back.

Barring that, sometimes a plane needs decommissioning, as it just ain't good to fly! Don't panic! As long as you stay subscribed to replacement tips, you have a lifetime warranty. Just drop us an email with a photo/video at hello@burstoralcare.com and our incredible flight attendants/customer service team will get you airborne again.



Contact support

Get on the intercom and chat with our divine customer service team at hello@burstoralcare.com

Random but important info

Our water flosser and compatible replacement tips are super tight. They're soulmates. No other tips but ours would do for the BURST Water Flosser!

When fueling up, only water or dental professional-recommended solutions will do. Don't take anything else off the beverage cart!

Lastly, remember that our water flosser is designed to clean gums only. No funny business! Use our product only for its intended purpose.



Water flosser ready for retirement? Give it the send-off it deserves – dispose of the battery properly with help from [burstoralcare.com/battery-disposal](https://www.burstoralcare.com/battery-disposal)

So what's next?

To infinity and beyond!

Ok, so we may not be literally going galactic, but we're sure going to aim for the stars when it comes to our new products. We will continue to expand our offering in order to provide total oral care solutions. We want your teeth to be treated like royalty. With the help of our free-wheeling, forward-thinking team of BURST Ambassador Dental Professionals, we've got some incredible pieces of kit coming to you soon.



May the fun
long continue...



**Feel like chatting with your
friendly neighborhood oral
care enthusiasts?**

Reach out to
hello@burstoralcare.com



.....
@f @BURSToralcare
#BURSToralcare
www.burstoralcare.com

BURST™ is a trademark
of BURST.USA Inc.

BURST.USA Inc.
340 S. Lemon Avenue #4658,
Walnut, CA, 91789